

# WTI Clinical Checklist

**Patient Name:** \_\_\_\_\_

**Laboratory Name:** \_\_\_\_\_

**Case Number:** \_\_\_\_\_ **Date:** \_\_\_\_\_

Checklist courtesy of: Dr. Stephen Wagner, D.D.S.

## Overview

- Wagner Try-In (WTI) is an evaluation of the final denture design.
- The WTI can be adjusted and equilibrated like any a conventional wax-up.
- Adjustments to the WTI should be made by warming wax and moving teeth or expressed in millimeters of change.

## Step-By-Step Checklist

- ☐ **Try in the WTI and adjust sharp spots or pressure areas.**
- ☐ **Does the WTI have good retention?** Yes/No **If NO**, take a wash impression of the WTI device for new scans using a light viscosity material.
- ☐ **Do the peripheral borders replicate the final impressions?** \_\_\_\_OK \_\_\_\_Needs adjustment  
If adjustments are needed, grind away over extensions.
- ☐ **Is the midline in position?** \_\_\_\_OK \_\_\_\_Needs adjustment  
If adjustment is needed, warm wax and move anterior teeth left or right to correct midline position.
- ☐ **Lip Support?** \_\_\_\_OK \_\_\_\_Needs adjustment  
**TIP:** Check by looking at upper teeth on wet/dry line of lower lip. Evaluate "SH" sound for lower teeth.  
If adjustment is needed, warm wax and move teeth facial or lingual to correct position.
- ☐ **Incisal Edge Position?** \_\_\_\_OK \_\_\_\_Needs adjustment  
**TIP:** Check whether teeth are too high or low, "F" and "V" sounds.  
If adjustments are needed, warm wax and tilt incisal of anterior teeth facial or lingual.
- ☐ **Is there a cant?** \_\_\_\_OK \_\_\_\_Needs adjustment  
If adjustments are needed, warm wax and move teeth left or right and up or down.
- ☐ **From a facial view, is the cervical of anterior teeth in proper placement?** \_\_\_\_OK \_\_\_\_Needs adjustment  
If adjustments are needed, warm wax and move teeth to lengthen or shorten their position.
- ☐ **If lengthening**, either scribe the desired position of the cervical of anterior teeth into the wax or provide a measurement from current position of cervical +/- in mm, relative to the High Smile Line. \_\_\_\_\_mm


☐ **Is there a preference for anterior tooth arrangement?**

**TIP:** Below are examples of different styles of tooth arrangements.

If adjustment are needed, warm wax and rotate individual teeth as necessary.

**Comments on arrangement:** \_\_\_\_\_

A- Rotate lateral incisors

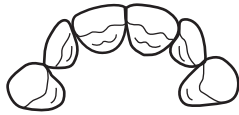

B- Flare central incisors



C- Small diastemata all anteriors



D- Move laterals to palatal



E- Rotate #9 with small overlap



F- Regular setup



☐ Once all adjustments are made, please have patient close in centric relationship, repeat multiple times to ensure accuracy of the VDO (Vertical Dimension of Occlusion). If the bite needs to be opened, squeeze the posterior bite blocks to increase vertical height. Again, have the patient close in centric, multiple times to ensure accuracy and take a bite registration.

**Indicate final Lucitone 199® gingival shade:**

\_\_\_ Original \_\_\_ Light Reddish Pink \_\_\_ Light \_\_\_ Dark Pink

**Indicate final tooth shade requested:**

\_\_\_ A-D, BW1, BW2, or Portrait® IPN® shade

**Check all of the following options desired:**

\_\_\_ Rugae \_\_\_ Post Palatal Seal \_\_\_ Exclude teeth

**ATTENTION CLINICIAN:** Provide precise notations for your laboratory and return this checklist with the WTI for additional adjustments and final design.